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Electric micro-needle instrument Instruction Manual

Basic facts of micro-needle 01

Nano technology
Nano wafer function: useful in makeup base foundation, skin moisturising and nourishing, suitable for daily maintenance, infiltration of various essence skin care products.
Low: base maintenance.
Medium: Zit, shallow pores, pregnant lips, infiltrate essence.
High: Mask, deep pox pores, foundation base.
Highest: Deep acne scars and deep pores.

Area of application
1.Skin care area:
Used in daily care, with the use of skin care products, to promote the efficient absorption of nutrients.
2.Skin repair area:
Used in acne, spots, scars, aging, etc., with the use of proper products to promote their rapid and efficient absorption to achieve better healing results.

Instruction for the use of electric micro-needles
1. Method of use with original base fluid or essence:
(1)Clean skin.
(2)Apply base solution / essence.
(3)Apply nano wafer vertically on the skin for 10-15 minutes.
(4)The excess essence left on the face can be gently pressed to allow it to absorb naturally.
2.Mask paste / eye patch:
(1)Clean skin.
(2)Apply mask / eye patch.
(3)Apply nano wafer to act vertically on the skin for 10-15 minutes
(4)The extra essence residue on the face can be gently pressed to let it absorb naturally.

Before Use - Frequently Asked Questions FAQs
1. How to clean.
Clean hands and face properly.Hand sanitizers or soap are usually used to clean the hands, and the face is cleaned with

Basic facts of micro-needle 02

facial cleanser or cleansing milk.
2. Whether the wafer should be sterilized.
The wafer is sterile and does not need to be sterilized.
3. The period for infiltration treatment.
It is normal to use it twice a week, and for sensitive skin should only be used once a week.
4. What kind of product is suitable for infiltration.
It is recommended to introduce disposable single usage products with ingredients of a certain effect or nature such as FGF, basic fluid (hyaluronic acid, VC etc), microcrystalline water light needles, and toning mask products.
***Before using strong effective ingredients, firstly do a test by applying a small sample behind the ear. If no allergic symptoms such as redness/swelling / itching then apply on the face please(stop using if allergic).**

During Use - Frequently Asked Questions FAQs
1.How to use the device.
All kinds of wafers are suitable for lifting and sliding, depending on the thickness and tolerance of the skin. Poor tolerance use point lifting, strong resistance and thick skin use sliding. Generally three areas namely eyes, forehead and mouth, and zygomatic or neck areas because of lesser subcutaneous fat, lifting is recommended. For the first time, it is recommended to use the lifting method. If there is no discomfort, the second time you can slide gently gradually improving your technique and strength control. For the initial few times it is necessary to find your touch. With the experience you can develop your own technique and usage of strength. on later treatment
2.How to control the timing and strength applied using the device.
For the initial usage control the device so that it is just below the wafer. Later you can apply tighter to the skin and develop skills through experience in handling the device and strength applied. Each time measuring the comfort of the skin as a gauge slowly applying more pressure and retracting if there is discomfort.

Basic facts of micro-needle 03


3. How to infiltrate various products
Essence:After applying the essence to the face, apply nano wafer vertically on the skin to activate infiltration. Can apply essence fully at a time or partially part by part.
Mask(non-washing mask):After the facial mask is used on the face,apply nano wafer vertically on the skin for 5-10 minutes.
4. Can you mix infiltration of several products at once?
Recommend that no more than two products should be infiltrated at a time, for sensitive skin should only infiltrate a single product at a time.

Frequently Asked Questions - After Use
1. Do you need to apply other skin care products after treatment?
It is recommended to use a moisturising or repair mask. If the skin is dry can be combined with moisturizing milk or cream.
2. Can touch water? Can you wash face or makeup the next day?
It is recommended to use the treatment in the evening and keep dry. If there are no symptoms such as red swelling and itching on the second day, normal skin care and makeup can be used. If you feel uncomfortable, please apply basic water supplement treatment.
3. Can go out in the sun, eat spicy seafood, alcohol?
It is recommended not to expose to the sun, eat less spicy food, avoid seafood, alcohol and other irritating foods within three days after use.
4. Can the nano wafer be reused?
The wafer cannot be reused. The wafer after expose will be oxidized and after use will be stained with remnants from the skin and follicle which will breed bacteria and it is therefore not recommended for reuse.
5. What if after treatment there is rash and itchiness?
It can be used daily for 3 consecutive days to heal moisturising mask and relieve sedation. in general it will generally improve significantly the next day. The fast and efficient absorption will achieve better repair effect.

Basic facts of micro-needle 04

How to apply FGF (Fibroblast growth Factor) freeze dried powder
Mix a bottle of lysase with FGF powder before use. The resultant FGF essence should be used within 3 days. FGF essence is more reactive and thus the skin can be susceptible to allergies. First smudge area with FGF essence wait for two days and if no allergic symptoms then apply infiltration treatment. If there is acne on the face, do not apply the treatment on the whole face. Use only in zit pores or scare areas. Do not use FGF with alcoholic products.
Advice on various problem skin use procedures and products usage
Dry, sensitive skin:
Before use, for the first three times introduced hydration products(hyaluronic acid solution moisturising mask) on the skin. After solving the hydration problems then introduce advance treatment.Test on a small area of the skin for allergies. If no redness, swelling or itchiness then apply on larger area or whole area. Apply non-washing moisturising products timely after treatment.
1. Spotted skin (freckles):
Suggest the use of Freckle essence in conjunction with water supplement products(such as whitening Freckle essence)
2. Acne problem skin:
Can be introduced with acne products, apply treatment carefully and within the treatment area. Avoid infecting other areas of the skin as acne is infectious.
3. Acne + Acne Case infiltration to be noted:
Can use acne products with FGF. First, use FGF on the acne area only avoid infecting other area of skin apply the treatment on the acne.
4. Specific instruction on pores treatment:
After cleansing, use the instrument to pick it up on the pore and then apply the FGF to the pit. Apply till redness and a slight pain will occur. Infiltrate twice a week. Daily apply FGF in the morning and night. Generally continued use in more than three months you can notice obvious improvement.

Description of Part 05

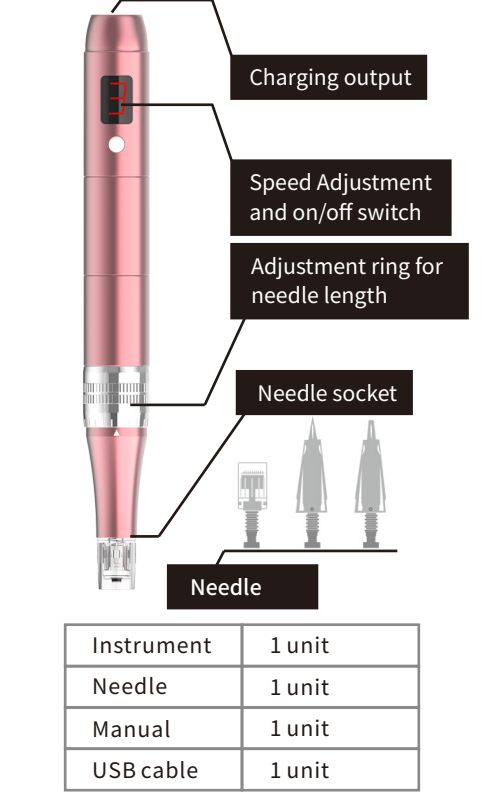Specification/Important facts 06

1. Name
Electric micro-needle
2. Functions
A micro-needle placed upon the skin between an extension and a retraction position to introduce a drug or cosmetic nutrients
3.Specifications
Product Model : WZ-103
battery : DC3.7V
Battery capacity : 350mA
Rated power : 1.5W
Operating voltage : DC3.7V
Charge voltage : DC5V
Net weight : 10g

Notes
Attention to note in usage of Electric Micro-needle:
01. DC5V charger must be used otherwise it may cause fire or malfunction
02. Do not use near high temperature, dust or humidity.
03. Proximity to the humidifier may cause malfunction.
04. Do not repair, disassemble or alter the instrument , it could cause malfunction.
05. Do not remove the power cord with wet hands, it could cause an electric shock.06. Use proper cable.
07. When you don't use it for a long time, pull it out of the electric socket.

Important facts 07

Skin Care Tips:
For proper skin care, it must choose products that are safe, effective, non-toxic, and suitable for one's own skin.
Things to note in usage:
1. Do not use if the disposable packaging is exposed and not sealed as it can cause skin infection.
2. Before treatment, the treatment area and hands should be clean.
3. When first used, it is recommended to rub skin care products on the inside of the wrist to perform a skin test. If there is itching and other allergic phenomena do not use.
4. Use moderate force in applying the wafer, do not intentionally press hard, and do not intentionally lift the wrist for contact.
5. About 10-20 minutes upon application, some will have a reddish skin phenomenon. This is normal and will disappear within 24 hours naturally.
6. Except for the wafer head, the rest of the instrument must avoid contact with fluids such as water so as not to damage the instrument.
7. The wafer is for single usage only. It is forbidden to share a wafer with another person.It is strictly forbidden to use the same wafer again otherwise it may cause adverse consequences such as skin allergies, redness, spots, and inflammation.
8. Discard used wafer and store the unused sealed wafer in well-ventilated and safe environment to avoid moisture and damage.Keep away from children.
9. During use, the wafer must be continuously in motion. Do not hover in the same area for 3-5 seconds exceeding twice.
The device is not suitable for the following groups of people:
People with acute fever or temperature 38 degree or more.
People with allergies or hemophilia.
People with melanosis.
People with heatiness.
People who often go to hospitals.
The pill man.
Under 12 years of age.

Method of Operation 09

1. Charge instrument upon receipt (DC5V)
2.Initial Preparation
(1) From the sanitary package remove the needle
(2) As shown in the diagram, attach the needle into the socket of the instrument till it touches the end.

(3) Turn clockwise to lock in needle.
(4) To remove the needle, turn anti-clockwise.
3. Adjustment of needle length
You can adjust the needle length by turning clockwise or anti clockwise. Clockwise to lengthen. Anti-clockwise to shorten.

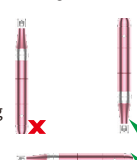
REGION	THIN SKIN	THICK SKIN
Forehead	0.25-0.5mm	0.5-0.75mm
Between Eyebrows	0.25-0.5mm	0.5-1.0mm
Nose	0.25mm	0.5mm
Eye Area	0.25mm	0.25-0.5mm
Cheek Bone	0.5mm	0.5-1.0mm
Upper Lip	0.25mm	0.25-0.75mm
Facial scars	1.0-1.25mm	1.5-2.0mm
Stretch Mark	1.5-2.0mm	1.5-2.0mm
Surgical & Burn scars	2.0-2.5mm	

4.Press and hold to turn on or off, short press to adjust the speed from low gear to medium or high

Operation Methods 10

Method of application
1. Cleanse skin:Daily makeup removal, face cleaning, can use toner for secondary cleaning.
2. Remove the wafer and install it on the top of the device.
3. Apply appropriate amount of basic fluid or essence to your face.
4. Activate device, vertical contact with skin surface, each contact point is recommended to hover 3 to 5 seconds only.
5. After the treatment is over, turn off the power and remove the wafer. Smear commonly used skin care products on the treated area. Recommended to use essence or cream. (Do a sterile repair mask before applying skin care products).
6. Wafer can not be reused. In order to avoid infection, please discard them after use.
Operational formula:
After the device is installed in place, it is directly in contact with the skin surface. It can visually detect a slight depression, and quickly and gently lifted up and down. According to different contact parts, adjust the strength.
Sliding method:
After installing wafer adjust device on the surface of the skin in contact vertically. A slight depression can be visually felt. Slide in a straight line and at the same time, make sure the wafer does not lose contact with skin surface. Do not hover or stay at the same area for more than 3-5 seconds and not more than twice.
Circle method:
Wafer installed and device held so wafer is perpendicular to the surface of the skin, You can visually sense the depression. Within a small circle move at a uniform speed making sure the wafer does not leave the skin surface. Again the same site is recommended to stay for 3-5 seconds and not more than 2 times.

Things to pay attention to 11

Usage of nano wafer
1.Cleanse face
2.Smear essence in need on the site to be treated.
3.Operation time is about 10-15 minutes. It is recommended that the first timer should apply gently and control time in about 10 minutes(touch contact to the surface do not apply too much force).During treatment, the same area does not exceed 3 seconds,Move back and forth about three or five times.
4.After treatment apply essence on the skin and wait for it complete absorption before applying a sterile repair mask.
Pre-use cosmetic product testing:
Before use, it is recommended to apply a small amount on the skin behind the ear. After 24 hours, if no red itch phenomenon then the product is applicable.
Post treatment:
1.After treatment, skin appears slightly red, slight prickle or pain, slight itch is a normal phenomenon. Can apply soothing cream or apply aseptic repair mask to relieve.
2.After treatment for 6-8 hours do not touch tap water (normally within 24-36 hours there should be no contact with water) as contact with bacteria makes it more susceptible to infection. You can wipe your face with saline water. Wafer treatment you can apply makeup after 24 hours while steel needle operation needs 3 days before light makeup.
3.After treatment, requires toning (morning / night application of aseptic facial/toning mask) and sun protection. Do not use skin products with alcoholic content and do not eat seafood within 3 days.
4.After two days you can do warm water cleaning and normal skin care.
Device Care
Be careful with the device and keep it from water or fluid. Use cotton stick to clean the device and needle preventing the remnants of any essence fluid getting into the device properly. Remnants can cause the device to jam.

Warranty/Maintenance 12

To provide you with better service, please fill the following form:

Warranty Period	One Year from the date of purchase		
	Yr	Mth	Day
Customer Information	Address		
	Name		
	Telephone		
Reasons for Repair			

We provide one year of warranty for the product from the date of purchase. If our product is faulty, we provide repair or replacement.
The warranty is void on the following clauses:
1.Damage due to failure to follow instructions provided by the operation manual.
2.Repairs done by unauthorized personnel.
3.Unauthorized tampering with the warranty card.
4.The serial number and model of product on warranty card does not match the product sent for repair.
5.Damage due to unforeseen circumstances.
Maintenance
1. General maintenanceMake sure the power is off and wipe the surface with dry cloth.Keep it away from direct sunlight or high temperature and humidity.
2. For equipment that is not used for a long time, make sure that it is safe to operate.