

MINI-RF

MINI RF Skin Beauty Device

Be sure to read this guide carefully
before you use this product



BEAUTY STAR

MINI RF Skin Beauty Device

USER'S MANUAL

Function and principles

When skin tissue's temperature rises over 40°C, blood circulation is expedited which promotes metabolism. In addition, hormone secretion becomes normalized due to the expedition of blood circulation. Cells activated with oxygen discharges toxic which prohibits secretion of melanic pigment. Using this mechanism, this device helps skin care as follows:

1. Managing wrinkles or skin rejuvenation

This device transfers high frequency RF energy to skin and stimulates the production of collagen in skin layer. In this case, stimulated collagen tissues are activated and thus helps prevent wrinkles due to stress and expedite skin rejuvenation.

2. Managing freckle and skin whitening effect

This device could alleviate freckles due to deteriorated blood circulation. High frequency RF energy transferred to skin increases blood circulation. This results in the reduction of melanic pigment secretion and thus whitens skin.

3. Managing pimple (Acne)

RF energy acting on the skin tissue generates heat which has an antibiotic effect. It normalizes sebaceous glands and regulates secretion of sebum which helps to reduce pimples (acne).

CONTENTS

Safety rules-----	2
Product specification-----	3
How does the Radio Frequency work-----	4
Massage Methods-----	5
Methods of Operation/Prohibition-----	6
Precautions -----	7
Function and principles-----	8

To our Customers

MINI RF Skin Beauty Device

Thanks for purchasing our high frequency technology skin care device. Please be sure to read all parts of this manual carefully for effective and safe use. You can also seek information about the principles of operation in our website. In addition, we provide updates for most recent information about this high frequency skin care device.



Precautions

Make sure to read the precautions before using the device and follow the instructions below:

- Use of gel to enhance stable effects for applying this device .
- Apply the gel(2-3mm thickness and 6cm by 6cm wide)on the area to be treated by around 2 times pumps. After treatment on gel-applied area, move to other areas and apply the gel as before, then restart the treatment.
- Please remove gel stains on the electrode of hand piece by soft napkin or cloth. Do not apply too much gel.
- Indicators on the device depends on device modes.
- Standby mode: when power is plugged into device, it automatically checks whether there is any error or malfunction. If there isn't any, it will be on standby mode in which output is minimized waiting for the next indication.
- Operation mode: when hand piece is contacted with skin at standby mode, clicking on the ON/OFF switch the device will start operating and RF energy will be emitted.



4 Pole RF Radio Frequency Head(Body)



6 Pole RF Radio Frequency Head(Body)



3 Pole RF Radio Frequency Head(Face)

Methods of Operation

1. Place this device on a flat floor for easy airflow.
2. Hand piece cable should be fully plugged into the hand piece.
3. Insert power cable into adapter terminal.
4. When the indicator is on, it means the machine is on standby.
When the indicator is flashing, it means the machine is in operation.
5. Apply 2-3 mm thickness of gels on area to be treated by around 2 times pumps.
6. After setting energy level. start operate referring to above diagram.
When gel is dried during operate it.
7. Move hand piece slowly over the area as shown in above diagram
but don't linger on a single area for too long.
8. After finishing treatment on gel applied area, proceed to treat other
areas as specified in the diagrams.
9. Remove gels on treated area after treatment and also gels remnant
on the hand piece and electrode
with cloth or napkin.
10. Applying ampoule or cream on skin after operation could enhance
skin treatment progress. It is recommended to apply UV block lotion
on skin after operation if you expose to the sun.

Prohibition

- Those who are equipped with cardiac pacemakers.
- Women who are pregnant and under breast feeding.
- Those who underwent face surgery..
- Those who had radiation treatment within the last 3 months.
- Those who have skin cancer around area to be treated. This device shouldn't be used at around the eyes, eyeballs, ears.
- Area with scar or wounds and skin with atopic or allergy.
- This product shouldn't be used for areas with skin problem or pustular acne.
- This product shouldn't be used for areas with itch and wound.
- This product shouldn't be used for those who are epileptic, mentally ill, pregnant, thrombosis patients, phlebites, and those who are equipped with pace maker, migraineur or patients with low or high blood pressure.

Specifications

Specifications:

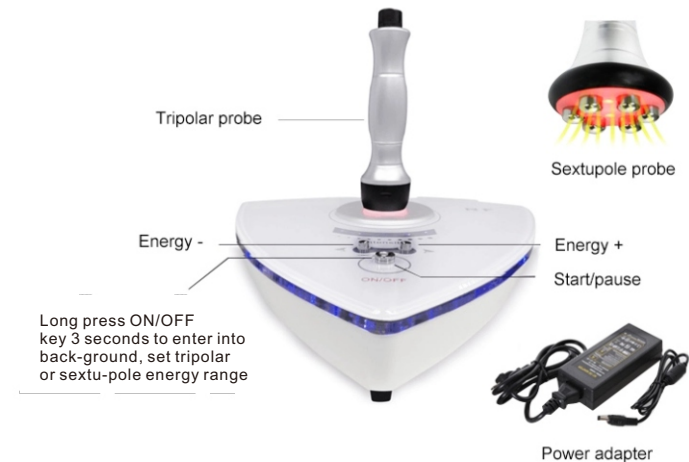
Energy: 1-15J/cm²

Power: ≤75W

Voltage: AC 110V-240V

RF Head: Three Head & Six Head

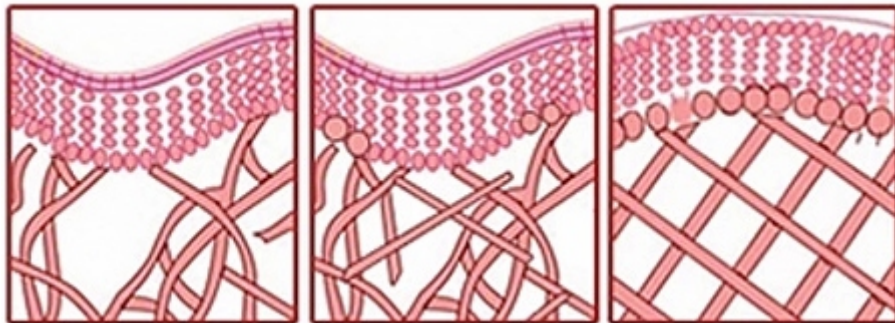
Main Device Size: Approx. 22 * 22 * 8.5cm






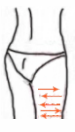






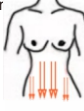
What is Radio Frequency and how does it work?

Radio Frequency (RF) is often used to reportedly tighten the skin. The (RF) energy penetrates the skin and stimulates the contraction of collagen and the production of new collagen tightens skin. Radio Frequency devices use a variety of delivery systems and penetrate to varying depths. The RF energy provided by our equipment, transferring a uniform volumetric heating effect into the deep dermis and underlying tissue. Heat is generated by the tissue's natural resistance to the flow of current within an electric field (impedance).

When the collagen and fiber in the derma layer heated to 45°C-60°C, the skin tightens at once and the collagen will constantly increase which can actually reach the effect of dispelling wrinkle forever and maintain youth.



Massage Methods

<p>Face (Treatment time: 30mins)</p>	<p>Follow the vertical direction to the ear shown by the hand gesture</p> 	<p>Follow the hand gesture clockwise and circle the face</p> 	
<p>THIGH (Treatment time: 30mins)</p>	<p>From down to upper pushing to the groin to dredge the lymph</p> 	<p>Horizontally along the thigh follow the hand gesture back and forth</p> 	<p>Also push by back and forth press motion to decompose fatness</p> 
<p>BACK (Treatment time: 20-30mins)</p>	<p>Push back and forth the bladder nerve 2 to 3 times.</p> 	<p>By lymph drainage system as illustrated, pull the toxin to the lymph node</p> 	
<p>ARM (Treatment time: 30mins)</p>	<p>From the inner elbow to the armpit follow lymph drainage system as illustrated</p> 	<p>From the elbow to the outer, do tightening gesture as illustrated</p> 	
<p>ABDOMEN (Treatment time: 20-30mins)</p>	<p>From small to large circles, in clockwise direction, circle around the navel</p> 	<p>Follow the lymph drainage system, as illustrated hand gesture, transfer the toxin to around the groin region</p> 	<p>At the belly region pull towards the groin region as illustrated from the two sides inwards</p> 