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MANUAL FOR HOME USE OF THERMAL
BODY&FACE LIFTING

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AN INSTRUCTION MANUAL

Contents

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01	Product Overview	09	BIPOLAR integrated head operation procedure
22	产品概述	30	BIPOLAR综合头操作步骤
02	Instrumental Advantages	13	Operating process of UNILARGE focus head
23	仪器优势	34	UNILARGE聚焦射频头操作步骤
03	Instrument Parameters	17	Use Contraindication
24	仪器参数	38	使用禁忌
04	Explanatory Figure for Product	18	Matters needing attention
25	产品图示	39	热拉提操作注意事项
05	System Interface	19	Instrument Maintenance
26	系统界面图示	40	仪器维护
		20	Packing List
		41	装箱清单
08	Operating Procedures	21	Warranty Card
29	操作步骤	42	保修卡


Product Overview

We are so appreciated you select our newest beauty apparatus for facial wrinkle reduction, and Thermal body&face lifting is the new generation of non-invasive skin firming and perfecting technology. Thermal body&face lifting is a non-invasive beauty apparatus for facelift that has passed China's Compliance Certification and CE Certification. Thermal body&face lifting Focus RF is the new generation of RF firming technology, and the accurate action on surface layer, superficial layer, middle layer and deep layer to achieve the purpose of wrinkles reducing, tightening and beautifying.

Before using the beauty instrument, please read this manual carefully and strictly adhere to the guidelines.

Instrumental Advantage

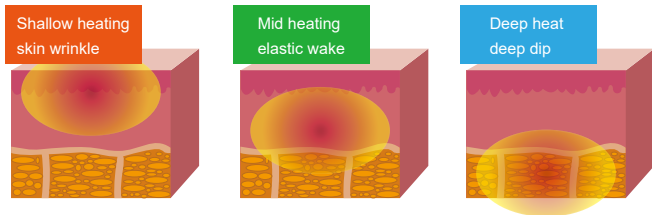"Thermal body&face lifting", a miracle instrument to make facelift without any surgical operation, which can effectively reduce the wrinkles, awaken the skin elasticity, lift the soft tissue structure, improve to firm the face and intensify the rejuvenation effect on the face, honored as "a spring of younger".

Instrument Parameters


Product name: Thermal body&face lifting		Model: H-020
Rate: 50-60HZ	Power voltage: 110~220V	Working voltage: 36V
Power: 80W	Instrument size: 25*25*10cm	Operating head: 2个


Explanatory Figure for Product


Explanatory Figure for System Interface

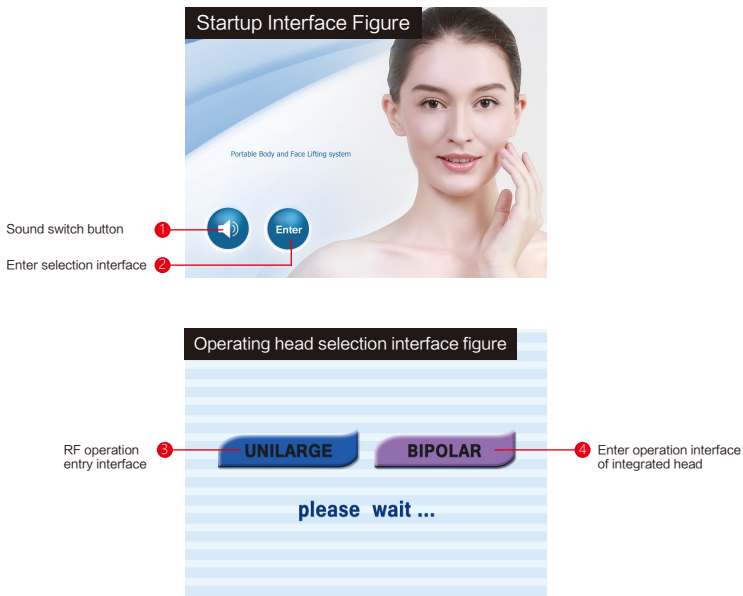
RF head operation interface figure

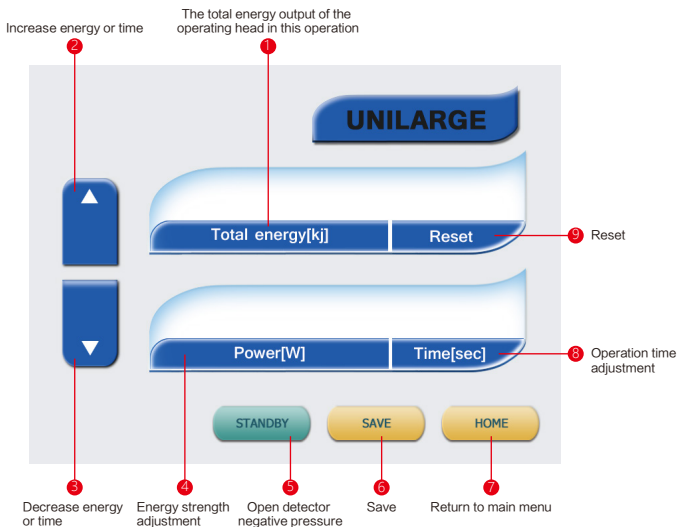
Integrated head operation interface figure

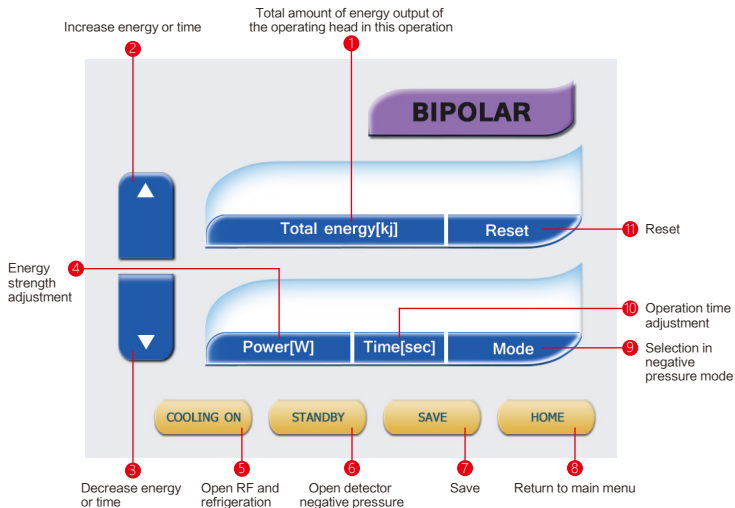
Operating Procedures

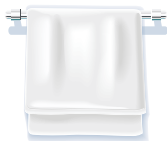
- Step1: Clean up the operation part firstly.

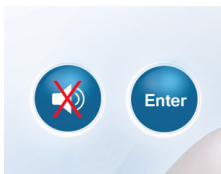- Step3: A small interface will occur on the display after opening the display (Figure 2), the speaker icon is the sound control. Click the speaker icon slightly, the speaker will display X, and X represents a closed sound of the instrument. Click ENTER to enter the main interface.

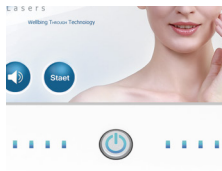
- Step2: Press the circular switch key at the top of the screen to open the display screen.

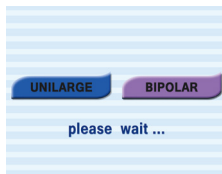- Step4: Two options are available after entering the main interface (interface figure). UNILARGE represents RF head operation interface. BIPOLAR represents to enter the integrated head (including three functions, i.e., negative pressure, refrigeration and radio frequency) operation interface.

BIPOLAR integrated head operation procedure

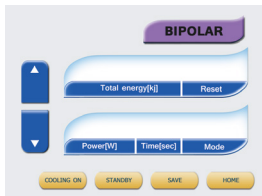Time[sec]

- Step5: We click BIPOLAR to enter the integrated head operation interface (as shown in figure).

- Step6: Then adjust the operation time, set the face operation time as 60min, and body operation time as 90 min.

Power[W]

1

Energy adjustment

2

Increase energy strength

COOLING ON

3

Open refrigeration and RF

STANDBY

4

Open the negative pressure function

Mode

5

Click item

- Step7: Click on the screen to select the energy adjustment, slightly click arrow to adjust the energy strength as 40, then click to open refrigeration and RF, click to open the negative pressure function. Finally click , four options are available for MODE, 1 represents that the operating head absorbs for 1s and releases for 1s, 2 represents that the operating head absorbs for 2s releases for 1 and 3 represents that the operating head absorbs for 3s and releases for 1s (in the operation process, the more time the operation head absorbs, the stronger feel it represents, vice versa). The mode should be selected based on the user's feeling. When the instrument is applied to the body, it is suggested to select the mode with the longest time for the body operation).

Power[W]

STANDBY

- Step8:Special medium is painted on the operation part.The customer can adjust the feeling in the operation process to adjust the energy strength, the energy will be increased gradually until the customer has the most comfortable feeling (it is strictly forbidden to operate the instrument from the highest temperature, or the skin may be injured).

Remark: The operation modes of integrated head and Focus RF are shown as follows:

The special medium we provided shall be used when operating the integrated head, the integrated head is provided with negative pressure suction, and the operation medium used must be of volatility, or it will block the air pump, resulting in abnormal use of instrument; in case of any operation failure caused by that the customer uses other medium, the customer shall bear the repair cost of the parts.

Neck

Lift and move the detector from the inferior position of neck to the inferior position of the chin along the neck, and move the chin to the lymph node position of the ears. Do remember that the detector shall not be operated on the laryngeal prominence.

Face

Move and lift (3-6 times) the detector to the position of zygomatic bone, and then to the position of ear lymph node from the corners of the mouth along a circular arc.


Canthus

Slowly lift the detector from the position of orbital bones to the temple, the energy should not be adjusted as too high, and it is suitable when the customer has a warm feeling, and move the instrument to the canthus along the orbital bones at a shape of semicircle. It is strictly forbidden to operate on the upper eyelid and eyeball.

Abdomen

Slowly move the detector to the navel from outside to inside centering on the navel, the horizontal operation should not be too large, and move the detector from the position of stomach to the cadera. (It is forbidden to operate it on the belly button)

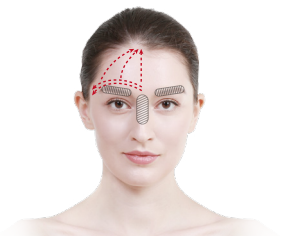Forehead

Lift the detector from eyebrows to the hairline, operate from the bottom up, then move and lift the detector horizontally to brow ridge around. It is strictly forbidden to operate the instrument on the eyebrows and bridge of the nose.


Leg

Move the detector from inner thigh to the outer thigh, from the knee to the leg end; one more operation should be conducted after the above cycle operation.

Arm

Move the detector from the arm downside to the outside, and slowly lift it to the outside (butterfly sleeve and bye-bye fat can be tightened)

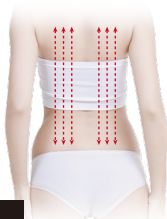
Back

Operate the detector from lower part of cervical vertebra to lumbar; repeat until cover the whole back. (Repeat 2-3 times for each action). Operate it on the back for scrapping detoxification.


- Step 9: Upon operation, you should clean the medium on the operating position, and apply pure hydrating mask and eye mask for 15 minutes.
- Step 10: Alcohol swab (subject to no drip) is used to wipe and sterilize the detector to clean medium residue on the detector and prolong the service life of the machine.
- For sensitive skin, it is recommended that you can use it every 15 days (integrated head mode 1 is used during operation, and it should not be too high for energy strength); for normal skin, you can use it every 5-7 days.
- **Attention:** The detoxification phenomena at different degrees will occur on the face after the operation: Such as: pimple, acne, acturc discascs, rubedo and itching are normal phenomena, which will be recovered to be normal on general skin within 1-3 days, and longer time for special cases. It is suggested to use hydrating mask every day for the above phenomena continually for one week. If the pimples or acnes occur after the operation, vitamin B6 ointment can be used for 3-5 times every day. As the machine can promote the skin blood circulation and accelerate detoxification of the skin, if the toxin is much in the skin, the pimples or acnes will occurs. If the skin is insufficient in massage care, lacking water or sensitive, rubedo, acturc discascs and itching will occur.

Attention: The bearing capacities of people's skin are not the same, thus, the energy strength used by the customer in the operation process is different, and the customer should adjust the energy strength of the instrument based on the actual feeling of your skin.

Operating process of UNILARGE focus head

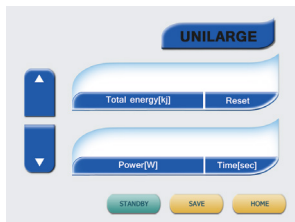- Step1:After completing the integrated head, click HOME to return to startup interface, click ENTER button to enter the main interface, then click UNILARGE to enter the operation interface of RF head.

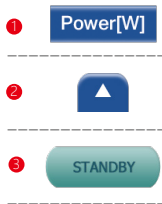
- Step2:Then adjust the operation time, set the face operation time as 60min, and body operation time as 90min.

- Step3:Click on the screen to select the energy adjustment, click arrow slightly, adjust the energy strength as 40 (the customer does not know whether it is suitable for the skin), click to enable the instrument to enter work preparation state.

Long press the button in operation

- Step4:Special medium is painted on the operation part.The customer can adjust the feeling in the operation process to adjust the energy strength, the energy will be increased gradually until the customer has the most comfortable feeling (it is strictly forbidden to operate the instrument from the highest temperature, or the skin may be injured).

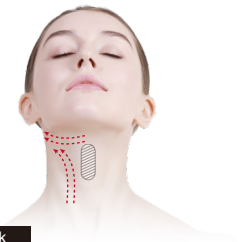Neck

Lift and move the detector from the inferior position of neck to the inferior position of the chin along the neck, and move the chin to the lymph node position of the ears. Do remember that the detector shall not be operated on the laryngeal prominence.

Face

Move and lift (3-6 times) the detector to the position of zygomatic bone, and then to the position of ear lymph node from the corners of the mouth along a circular arc.


Canthus

Slowly lift the detector from the position of orbital bones to the temple, the energy should not be adjusted as too high, and it is suitable when the customer has a warm feeling, and move the instrument to the canthus along the orbital bones at a shape of semicircle. It is strictly forbidden to operate on the upper eyelid and eyeball.

Abdomen

Slowly move the detector to the navel from outside to inside centering on the navel, the horizontal operation should not be too large, and move the detector from the position of stomach to the cadera. (It is forbidden to operate it on the belly button)

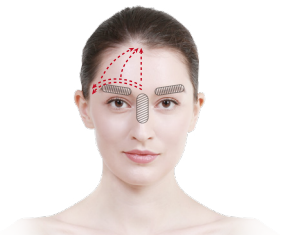Forehead

Lift the detector from eyebrows to the hairline, operate from the bottom up, then move and lift the detector horizontally to brow ridge around. It is strictly forbidden to operate the instrument on the eyebrows and bridge of the nose.


Leg

Move the detector from inner thigh to the outer thigh, from the knee to the leg end; one more operation should be conducted after the above cycle operation.

Arm

Move the detector from the arm downside to the outside, and slowly lift it to the outside (butterfly sleeve and bye-bye fat can be tightened)

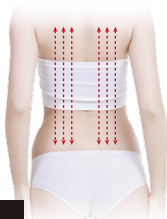Back

Operate the detector from lower part of cervical vertebra to lumbar; repeat until cover the whole back. (Repeat 2-3 times for each action). Operate it on the back for scrapping detoxification.

- Step 6: Upon operation, you should clean the medium on the operating position, and apply pure hydrating mask and eye mask for 15 minutes.
- Step 7: Alcohol swab (subject to no drip) is used to wipe and sterilize the detector to clean medium residue on the detector and prolong the service life of the machine.
- For sensitive skin, it is recommended that you can use it every 15 days (integrated head mode 3 is used during operation, and it should not be too high for RF energy); for normal skin, you can use it every 5-7 days.
- **Attention:** The detoxification phenomena at different degrees will occur on the face after the operation: Such as: pimple, acne, acturc discascs, rubedo and itching are normal phenomena, which will be recovered to be normal on general skin within 1-3 days, and longer time for special cases. It is suggested to use hydrating mask every day for the above phenomena continually for one week. If the pimples or acnes occur after the operation, vitamin B6 ointment can be used for 3-5 times every day. The machine is effective to promote blood circulation of skin and boost the toxin emission at the bottom of skin. Heavy toxins remaining under cover of skin will spur the occurrence of pimples or acne; If the skin is insufficient in massage care, lacking water or sensitive, rubedo, acturc discascs and itching will occur.

Attention: The bearing capacities of people's skin are not the same, thus, the energy strength used by the customer in the operation process is different, and the customer should adjust the energy strength of the instrument based on the actual feeling of your skin.

Use Contraindication

People with the following symptoms shall be cautious when using the instrument and ask for advices from the doctor or professional; the details are as follows:

- In case the operation part has been filled with artificial fillers, such filled part shall be avoided when operating the instrument.
- The customer should tell to doctor whether she/he has injected hyaluronic acid or collagen within half year;
- The customer shall not wear any metal articles in the operation;
- People with heart diseases or wearing the pace maker are not suggested to conduct any operation.
- The customer is not suggested to make up on the day accepting the operation, or the face should be cleaned up before the operation;
- The operation shall not be operated if the female is in the period;
- In the operation, the detector shall not be operated on the position with wound or purulent acne, within orbital bones or Adam's apple;
- This instrument is not suitable for the patient allergic to metal.
- Although sufficient evidence is not available to certify that this instrument can influence the patients with tumor and active skin disease, it is not suggested that such patient to use this instrument.
- Although direct evidence is not available to certify that this instrument has any effect on the pregnant, it is not suggested the pregnant to use this instrument.
- Although it is not certified that the ultrasonic wave may result in osteoporosis of the rib, it is not suggested to operate the instrument on the rib.
- If the face is in an allergic period, this instrument shall not be operated within one week after the completion of the allergic period.

- The lowest level of energy shall be adopted when using this beautifying instrument, and the energy can be adjusted and increased gradually based on their actual feelings, until they have the most comfortable feeling. It is forbidden to open the energy to the largest level, the bearing energy of each person is not the same, thus the energy strength will be different in specific operation, and if the instrument is started with the highest energy level, it may injure the skin.
- People not suitable to use the instrument: 1. Those do not carry out any care and maintenance; 2. Those have a sensitive and extremely dry skin; 3. Those have pustules and acne. After the operation, detoxification phenomena may occur on the face, such as pimple, acne, acturc discascs and rubedo, which are normal. For such conditions, maintenance and hydrating mask will be effective.

Note: This instrument should be matched with our dedicated medium instead of other gel or essential oil when using; any problem incurred due to using other mediums will be considered as individual responsibility.

Notes to Operate Thermal body&face lifting

- Hot water shall not be used to wash face within 3 days; (warm or cold water is available)
- The moisturizing and sun protection of the operated position shall be strengthened;
- Hot spring shall not be used within 3 days;
- Make-up shall not be conducted within one day after the operation, light make-up is available after one day, and heavy make-up shall not be conducted after 4 days;
- It is suggested to use 4 masks at least, and drink water not less than 2000cc (8 cups of water).
- Spicy food shall be avoided as much as possible, smoking, drinking and staying up shall be forbidden, fruit and vegetables shall be eaten as much as possible, and the meals shall be light.
- It is not suggested to carry out strenuous exercises at the current day using the instrument.

Equipment Maintenance

- Please clean the instrument with clean water and stored properly after each operation.
- Plug with pin is a must before using the instrument, and it must ensure that the power socket of the instrument is really grounded.
- Ensure the instrument voltage is suitable. We suggest the user to increase the voltage-stable power matching with the power if the local power voltage is not stable.
- For treatment effect of the instrument and ensure the normal service life, please use the appointed accessories provided by the originating factory.
- The instrument can neither be placed in the damp place or around water source; nor be exposed to direct sunshine.
- The instrument cannot be close to the strong heat source, if so, it shall affect the service life and normal use of the instrument.
- Current loop guard shall be used in winter or dry weather, the position contact the skin shall be provided with a certain moisture to avoid the poor contact, and wiping with water-based product or water can reach the necessary effect.
- The medium, such as water, oil, etc. cannot leak in the instrument, and the instrument cannot be knocked, thrown or heavy falls in case of damage the instrument.
- Please clean the detector timely after the operation and cut off the power supply.
- The warranty period is one year. The instrument is maintained and replaces components for free under the precondition of normal use during the warranty period.
- We provide lifetime maintenance and technical support.
- The instrument shall not be cleaned and wiped with any articles containing ethyl alcohol.

Packing list

1

Host×1

2

RF
operating head×1

3

Integrated
operating head ×1

4

Power wire×1

5

CD×1

6

Certificate of approval×1

7

Instructions×1

Warranty Card

Maintenance Record Sheet

Product name: Thermal body&face lifting Purchasing date:

One year of warranty period (since the purchase date)

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