


D-ICE CARE HAIR REMOVAL USER MANUAL

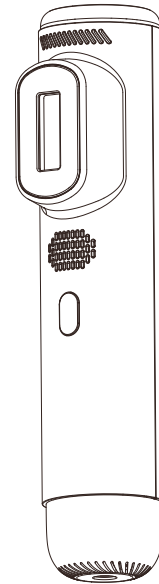

D-ICE CARE HAIR REMOVAL

Table of Contents

Product description-----	02
Safety precautions-----	03
Production features-----	05
Introduction of display-----	06
Description of Work Mode Switch-----	06
Ice-cool Care Mode-----	07
Skin sensitivity test-----	08
Operation steps-----	08
Course of treatment-----	13
Operation Method of Skin rejuvenation function-----	14
Suggestion on Operation Interval of Functions-----	14
Adverse Reactions and Treatment-----	15
Packing list-----	15
Specification Parameters-----	16
Troubleshooting-----	17
Applicable Part-----	18
Frequently Asked Questions-----	19
Operation record-----	22

Read this Manual and understand all instructions and precautions before operating the product

Production description

The ICE-Cooling care hair removal device (hereinafter referred to as the hair removal device) is a beauty product featuring hair removal, skin rejuvenation and ice-cold compress. With Home pulsed light technology, the hair removal device integrates many optical technologies applicable to the field of beauty and health, which is one of the professional and practical products in this field.

Based on the ice-cool technology, the hair removal device with ice-cold compress functions. When the ice-cool care mode is enabled, it can reduce the excessive heat generated on the skin by the photon irradiation and do ice-cold compresses during hair removal. It can effectively prevent burns and aches on the skin, thus giving appropriate care to your skin. Customer's comfort and safety always come first in the design of the hair removal device, which is manufactured in line with high quality standards.

Applicable customer base

The hair removal device is an over-the-counter instrument used to remove excess hairs for adult men or women. The treatment is a totally physical therapy.

Suitable skin color

The use for black skin is prohibited. Other skins are applicable. It is recommended to remove the hair on brown skin at the lowest level, which can be adjusted to the corresponding level after your skin is adaptive to the pain.

Safety precautions

Contraindications

- 1.Children and elders;
- 2.People with systemic inflammation;
- 3.Photosensitive skin and people using photosensitive drugs
- 4.People using a life support unit such as an implanted medical electronic machine (e.g. a heart pacemaker) and an artificial heart -lung. The product may result in incorrect operations of the electronic medical equipment.

Safety Warnings

- 1.It is forbidden to disassemble the hair removal device without permission given the high voltage inside the product.
- 2.Never cover the ventilation port when using the product, which is used for heat radiation and temperature reduction for the product.
- 3.The product shall be kept away from moisture and water.
4. Do not use the product in case of damages to plug or the power cord.
- 5.It is forbidden to operate the instrument in a flammable and explosive environment.

Precautions


- 1.The product is not recommended for women while pregnant or trying to get pregnant or nursing since this product have not tested among these women group.
- 2.The instrument shall not be used continuously for more than 25 minutes. If the treatment area is relatively large, it can be divided into several times. (Be sure to power off the product and unplug the power after 2 minutes of use. Continue to use after 10 minutes.)
- 3.The use of the product is not recommended after a sunbath, when a sunbed is used, or when a sunless tanner is used.

4. Do not use the instrument for positions with tattoo, or positions with dark gray or black spots. Otherwise, the skin may be uncomfortable due to absorption of a large amount of light energy.
5. Its forbidden to use the product for the ears, neck, nipples, genitalia or positions around the anus, because these positions are sensitive and susceptible to injury.
6. Irradiating during each treatment at the same area shall not be more than 3 times to prevent causing uncomfortableness to the skin due to excessive heat.
7. Suspend the use of the product if the skin becomes red and swollen and the condition persists for a minimum of 24 hours. In this case, contact a doctor or consult our customer service personnel immediately


Specification for Safety Operations

1. It is forbidden to aim the lamp at vulnerable parts such as the eye. Wear goggles whenever possible. Do not observe the light or reflected light from the light outlet of the lamp even if you are wearing goggles.
2. Make sure that all metal materials, such as finger rings and bracelets, are removed from you body before using the product. Otherwise, the remover may fail.
3. Make sure that all metal materials, such as finger rings and bracelets, are removed from you body before using the product. Otherwise, it may cause equipment failure.
4. Always keep away from children during operations. The product shall be properly kept at locations not accessible to children.

Production features

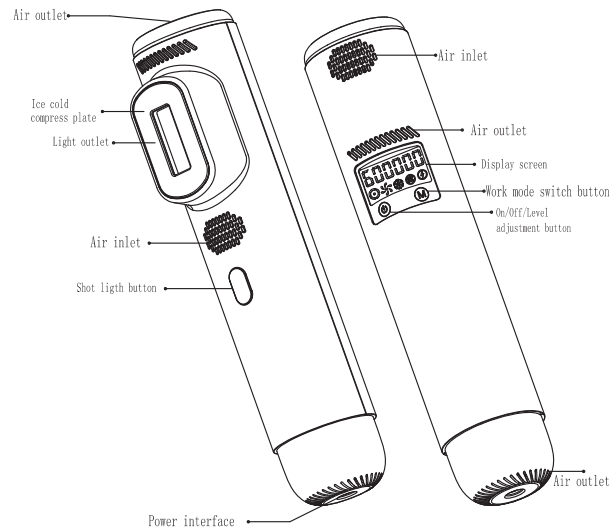

Introduction of display

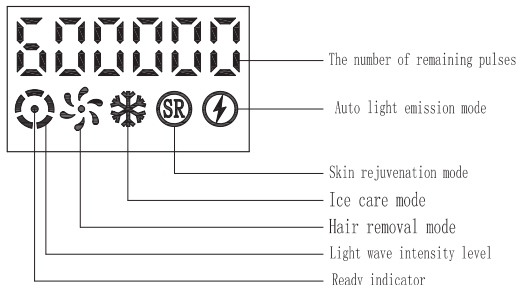

Description of Work Mode Switch

1. Refer to the "Operation steps" below to turn the device on.
2. Switching between the hair removal and skin rejuvenation modes: Press the "Work mode switch button"; SR indicates the skin rejuvenation mode, otherwise is the hair removal mode


600000	Beginning pulses is 600000, when the battery power getting low, it will be flash
•	Ready indicator, when it flash, means can start shot light
⦿	5 levels of light density
🌀	Normal hair removal mode (default mode)
🌀❄️	Ice care & hair removal mode
🌀❄️Ⓜ️	Ice care & skin rejuvenation mode
⚡	Automatic light emission mode

⚠️ Precautions

1. After the ice-cool care mode has been working for a period of time .it cannot be directly changed into the non-ice-cool mode to prevent burning your skin. If you want to change it to the non-ice-cool mode, you must turn off the device and wait for 20 minutes to dissipate heat, and then turn on the device to adjust it to the non-ice-cool mode.
2. After the ice-cool care mode is turned off, the heat will be transferred to the ice-cold compress panel. At this time, keep your skin away from the ice-cold compress panel. Turn the product off, and wait for 20 minutes to dissipate heat.

Ice-cool Care Mode

When the hair removal device is in the manual light emission mode (the default mode after starting up). Under the following two work modes: 🌀❄️, 🌀❄️Ⓜ️, you can carry out ice-cold compress care without pressing the shot light button.

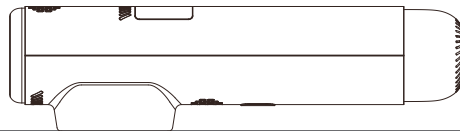

2. Attach the ice-cold compress panel vertically and closely against your skin, and start the ice-cold compress care.
3. Tips: In the process of ice-cold compress care, it is necessary to press the "Level adjustment button" once within 5 minutes to keep the hair removal device running. Otherwise, it will automatically shut down. In this case, just hold the "On button" down for 2 seconds to start it up

Quick start

Before operating the hair removal device for the first time, select the energy level suitable for your skin and the body part where you want to treat. The hair removal device has 5 energy levels from Level 1 (lowest) to Level 5 (highest). Level 1 is for skin rejuvenation, and Level 2-Level 5 are for hair removal

Skin sensitivity test

If you are a new user of this product, or if you are using it on some body part for the first time, first select a small area and perform the skin sensitivity test as follows:

1. Turn on the product with reference to the "Operation steps" below.
2. Select Level 1 or Level 2.
3. Operate the device on the selected skin area, following the instructions and the precautions above.
4. Check the reaction of your skin after 30 minutes of skin test.
5. If the skin tested is normal, you can continue to use the device.
If the skin appears red and swollen, carry out the skin test again 24 hours later.

Operation steps


1. Before use, please shave all hairs around the treatment area, clear up and wipe up first to ensure that hair follicle underneath the skin can absorb energy to the maximum extent possible, so as to achieve the best treatment effect.




Tips

It is recommended to shave the hair in the following two methods:

- 1). Lubricate the area to be shaved with shaving cream or mild shower cream
 - 2). Alternatively, cover the area to be shaved with a hot towel for a few minutes, and gently shave the hair immediately, which will make it easier to remove the hair.
2. Put the plug at the adaptor tail into the power interface of the hair removal device

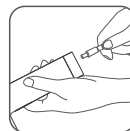


3. Connect the adaptor to the power outlet; the numbers on the display screen will light up, after about 3 seconds will turn off if no operation

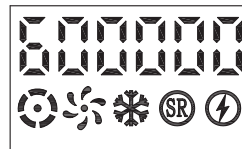

4. Press and hold the "On button" for 2 seconds to start the hair removal device. By default, it will enable the hair removal and Ice-cool care mode. Check the working conditions with reference to the "Introduction of display" and "Description of work mode switch" above.
5. This device is designed with 10 levels from Level 1 to Level 5. Level 1 is enabled by default upon its start-up. Press the "Level adjustment button" to adjust the level. Start operation from the low level. If you feel very comfortable, gradually turn it to the corresponding level

suitable for your skin

6. Put on the goggles (or sunglasses) supplied with the device to prevent injury to your eyes.

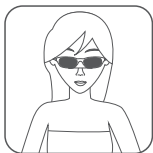

7. Setting of manual light emission and auto light emission

Manual light emission mode The default mode is manual light emission mode after the hair removal device is turned on. Press the "Shot light button" one time, light is emitted one time.

Auto light emission mode (for hair removal in large areas)

After the hair removal device is turned on, press and hold the "shot light button" for 2 seconds; the ⚡ on the display screen will light up, and the hair removal device will switch to auto light emission mode. You only need to press the lamp closely against the skin, and the hair removal device will automatically emit light. Attention: When the hair removal device is in the Auto light emission mode, it is strictly forbidden to touch the lamp with body parts that do not need to be treated. It is strictly forbidden for others to touch the lamp.

Switching the manual light emission mode

In the auto light emission mode, press the "shot light button", the ⚡ on the display screen will turn off, and the hair removal device is switched to manual light emission mode.

8. Put the lamp close to the position to be treated, and try to hold the lamp vertical to the skin:

- ①. When the hair removal device is in the manual light emission


mode, the readiness indicators flashes, press down the shot light button to irradiate. Let the lamp to cover the area needed to be treated to complete hair removal. Move the lamp to the next area needed to be treated, and close it to the skin. When the readiness indicators flashes, press down the shot light button to irradiate. Repeat these steps repeatedly. It is suggested to irradiate each area 2-3 times, not more than 3 times.

②. When the hair removal device is in the auto light emission mode, light will be emitted without pressing the shot light button. The lamp covers the area needed to be treated to complete hair removal. Move the lamp to the next area needed to be treated, and close it to the skin. The hair removal device emits light. It is suggested to irradiate each area 2-3 times, not more than 3 times.

9. After operation is completed, long press the Off Button for 2 seconds to turn off the hair removal device. When the hair removal device is in the ON state, it will automatically turn off if no operation is carried out within 5 minutes.
10. Unplug the adapter, and remove the adapter from the power outlet.
11. Wait for 20 minutes until the surface of the lamp is cooled; wipe the light outlet of lamp with a non-alcoholic wipe or a soft tissue.
12. After cleaning is completed, put the hair removal device and adaptor back to the packing box. Put the hair removal device in a dry and well-ventilated place, which should be kept out of reach of children.
13. Pamper the treated skin
- After 20 minutes since completion of operation clean the treated skin, and you can use the moisturizing skin-care products to moisturize and repair the treated skin. If you feel uncomfortable, apply cold to the skin


14. Charge the battery:

Plug the adaptor DC pin into the charge socket as show bellow
When charging, the charge indicator flash;
when full charge, the charge indicator light

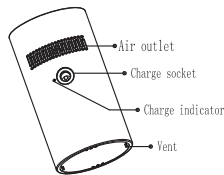

15. Using battery:

Insert the battery into the main unit, and the can operation as normal.
When the battery power getting low, the remaining numbers on the display screen will flash, please charge it on time

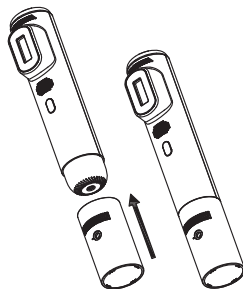

Attention:

- Apply the aloe vera gel or similar moisturizer to the area after treatment. Do not apply the irritating cosmetics immediately.
- Before operation, ensure that the area to be treated is clean. It is recommended to clean the skin with cold water. The treated skin cannot touch water within 2 hours after operation.
- Do not have a sun bath, warm water bath, steam bath, etc. within 48 hours after operation.
- Try to avoid exposing the treatment skin to the bright sun within 2 hours after the treatment. Good sun protection measures before outdoor activities. It is best to treat at night.
- Do not swim in the disinfected swimming pool within 48 hours after the treatment

Course of treatment

First month	Use it once every 3-5 days	The appearance of hairs or skin will not change obviously. Please be patient.
Second month	Use it once every week	The hair begins to become thinner and softer, and when taking bath, you can find that the black hair follicles will fall off, which is the happy surprise.
Third month	Use it once every 15 days	You can see that new hairs are basically suppressed, and will not grow again. New hairs can be pulled out gently. Only a few hairs grow out.
Fourth month	Use it once every 30 days	There is almost no hair; you complete long-lasting hair removal and get perfect skin look.

Remarks

The hair removal device is used to remove hair permanently according to period arrangement. Some treatment parts do not need follow-up maintenance. Some treatment parts need to be maintained for different times such as 6 months, 9 months and 12 months.

Operation Method of Skin rejuvenation function

- 1).After start-up, short press the Work Mode Switch Button to adjust it to the skin rejuvenation work mode.
- 2).Press the level adjustment button to adjust it to level 1, and then align the lamp with the area requiring skin rejuvenation for irradiation. It is suggested to irradiate each area for one time, and carry out operation 1-2 times every week. Insist on operation for 2-3 months, you can see obvious effect.

Suggestion on Operation Interval of Functions

Each position can only receive photon irradiation for one function in one day. If the same position requires irradiation for the other function, operation should be conducted on the next day.

Adverse Reactions and Treatment

All beauty treatment ways may have a certain risk. Similarly, if the hair removal device is not used properly, it may cause side effects. Side effects that may occur and suggested treatment ways are as follows:

Pain	Each person has different sense of pain due to different constitution. It is rare to feel uncomfortable. If you have slight burning sensation, which may last 1 minute at the most. Most adults and older children can bear this kind of discomfort.
Minor redness or slight burning sensation	It is normal. It is because that hair removal has not been done or hairs have not removed completely before using the cosmetic instrument.
Hair burnt smell	It is normal. It is because that hair removal has not been done or hairs have not removed completely before using the hair remover.

Packing List

Main unit	1	Adaptor	1
Battery (option item)	1	Goggles	1
Shaver	1	User manual	1


Specification Parameters

Dimension (main unit)	222x48x61mm
Weight (main unit)	270g
Irradiation area	3cm ²
Power supply	Input: 100-240V~50/60Hz Output: DC12V 4A
Battery	DC7.4V 2000mAh
Rated power	45W (max.)
Operating environment	Temperature: 5-40°C Relative humidity: 20-90%, without condensation
Storage and transportation environment	Temperature: -20 - 55°C Relative humidity: 5-95%, without condensation
Guarantee period	One year (for main unit only)

Trouble Shooting

Syndrome	Description	Solution
No display	The power cable is not properly connected.	Connect the power cable.
Startup failure	The On button has not been held for more than 2 seconds.	When starting, be sure to hold the On button for at least 2 seconds.
The light is not emitted when the shot light button is pressed	1.The lamp does not touch the skin. 2.Electric power storage has not been completed. The readiness indicators is not flashing.	1.Re-touch the skin. 2.Wait until the readiness indicators has flashed, and then press the button.
Automatic shutdown	The device is shut down automatically if no light is emitted within 5 minutes	Hold the On button for 2 seconds to start up the device.
The remaining number flash on screen	Battery low power	1).Use adaptor 2).Charge the battery
ERROR display on screen	System error	Cut off the power, cool down the main unit fully

Applicable Part Note: Do not get too close to your eyes when you are using the device!

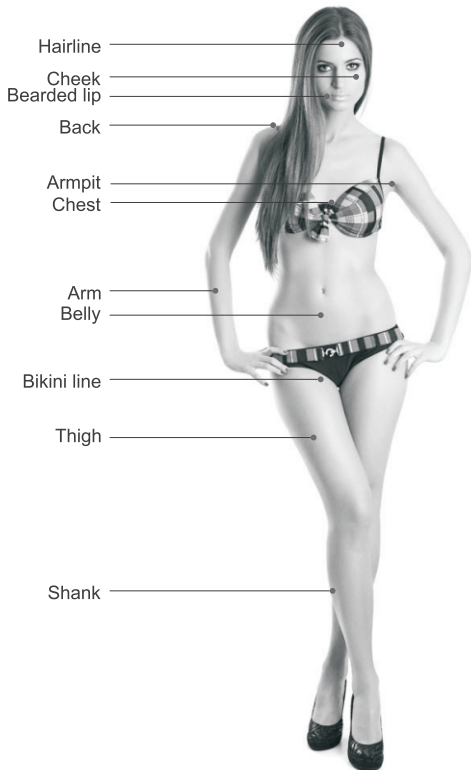

Frequently Asked Questions

Q: Which is the hair removal mode? Which is the skin rejuvenation mode?

A: Press the "Work Mode Switch" Button, it is the skin rejuvenation mode when the display screen shows SR. It is the hair removal mode when the display screen does show HR.

Q: Does level adjustment apply to ice-cool care effect?

A: Level adjustment is only for energy intensity, and has no effect on ice-cool care.

Q: Is dehairing required before each operation?

A: Not necessary. When you use the device for the first time, the thicker hair should be shaved once. When the hair has grown to approx 5 mm, and the use effect will be affected, the hair should be re-shaved one time.

Q: Why are not my hairs removed immediately after irradiation, or are not my hairs changed after irradiation?

A: It is a permanent hair remover. Hairs cannot be removed permanently until you have used the device for one period. It can only repress hair growing and slow down hair growing speed at the beginning. The hair enter the sleep state under continuous action. Thus, permanent hair removal effect is achieved.

Q: I have used the device for a few times. Why is the effect on my skin not as good as that of others, or has no effect achieved?

A: Everyone's constitution is different. Hair varies in thickness and color. The device is used on different positions. Therefore, use effect periods are different naturally. The use period is an average value. Detailed period may be different slightly. If no effect has been achieved after you has used the device for more than one

month, be sure to contact us. Absolutely, you has not operated the device properly.

Q: Which level is more suitable?

A: No uniform standard is used for this. It depends on each person's skin endurance (the darker the skin is, the lower level should be used). You should begin to use the device from lower level, and then use higher level gradually. The higher the level is, the energy is greater, the use effect will be more obvious, and the hair removal period is shortened correspondingly. In contrary, the lower level is used, the hair removal period is delayed correspondingly. It is suggested to carry out test irradiation on a small part of your skin for 1-3 times before hair removal on each position.

Q: Why does a hair burnt smell occur after I use the device?

A: Hairs have not been shaved completed during use. The melanin will interact with light wave to generate heat. Therefore, hair burnt smell will occur. It is normal.

Q: Will it hurt?

A: During use of hair removal device, light wave will generate heat when it interact with melanin. Therefore, some slight pain will occur. Everyone's skin sensitivity is different, pain tolerance is different. Therefore, everyone's feeling is different. If you feel pain, you can begin from the lower level and then use higher level gradually to let your skin adapt to it.

Q: How to shave hairs? Which hairs should be shaved off?

A: Armpit hair and hairs at the positions such as bikini line leg should be shaved before use of the device (It is suggested to shave off the hair with a length of more than 1 cm). Before shaving hair, it is suggested to use the following two methods to remove hair more easily:

①. Lubricate the position to be dehaired with shaving sream or gentle shower cleanser;

②. Or use a hot towel to cover the position to be dehaired for a few minutes, and then gently shave the hair. In this way, it is more easily to remove hairs.

Q: I want to remove hair more quickly. Can I use the device everyday?

A: It is not suggested to use the device everyday. Hair has a growth cycle. You just need to use the device based on the normal use cycle. You can use the device one time every 3-5 days in the first month. After the hair growing speed is slow down and the effect is obvious, you can gradually reduce use times.

Q: Why does my skin turn red and little red spots occur after I have used the device?

A: It is because that the skin is hurt by shaving too much based on many times of tests. It is suggested to use the following two methods to remove hair more easily:

Lubricate the position to be dehaired with shaving sream or gentle shower cleanser;

Or use a hot towel to cover the position to be dehaired for a few minutes, and then gently shave the hair. In this way, your skin will not be hurt. It is best to use the electric shaver to shave hair if possible. Everyone can test it on a position without shaving to see if this kind of situation will be resulted in.

Q: Why I cannot see effect after I have used the device several times?

A: The positions used with the device are different and hair thickness is different. Therefore, the effect will be different. In general, armpit hair and men's beard are more difficult to be removed due to male sex hormone and gene, the use effect will occur more slowly. The use effect is also related to the use level. The permanent hair removal effect will be achieved as long as you insist on using the device.

Operation record

Function selection: Hair removal Skin rejuvenation

Operation site _____ Level _____

Number	Date	Remark
1		
2		
3		
4		
5		
6		
7		
8		
9		
10		

Operation record

Function selection: Hair removal Skin rejuvenation

Operation site _____ Level _____

Number	Date	Remark
1		
2		
3		
4		
5		
6		
7		
8		
9		
10		

