

MANUAL

SKIN CLINIC



SKIN CLINIC

SKIN CLINIC 6 steps -----	02
PRODUCT PARTS -----	03
DEVICE FUNCTIONS -----	04-05
PRODUCTION DIRECTION -----	06-08
BATTERY/CHARGE TIME -----	09
STORAGE&CAUTION -----	10

COMPONENTS

○ CHARGER



○ DIRECTIONS



6 STEPS



STEPS 01 VIBRATION

Micro vibration massage for enhanced effect.



STEPS 02 HEATING MASSAGE

Heat massage helps to open pores to enhance function.



STEPS 03 NUTRITION

Supports cosmetics to penetrate deep into skin.



STEPS 04 FACIAL CARE

+,- current flows to stimulate and treat skin.



STEPS 05 BODY CARE

+,- current flows to stimulate and treat body.



STEPS 06 LED LIGHT CARE


LED care massage to enhance and support skin.

i PRODUCT PARTS

SKIN MANAGER


01 SKINIRON Device

o PRODUCT

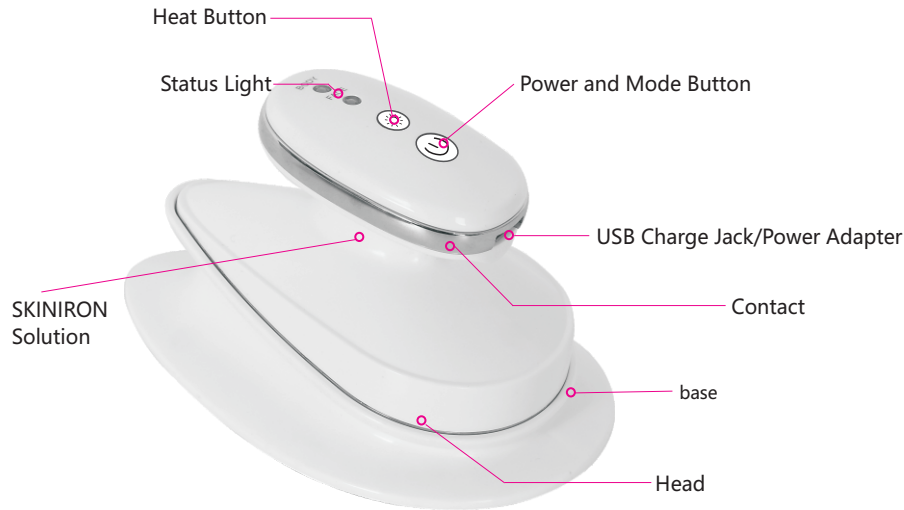

SKINIRON SOLUTION SPECIFICATION		
Rating: DC3.6V 2.600mA	Consumption:8W	Battery Life:approx 80mins
Size:L115 X H65 X W65(mm)	Weight:130g	Battery Charge Time:approx 3hrs

o PRODUCT STAND

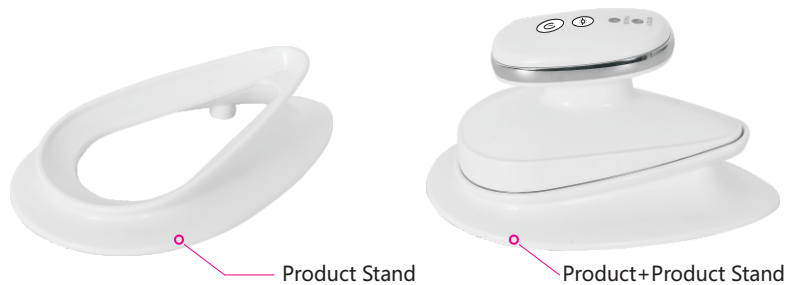

o DEVICE FUNCTIONS

SKIN MANAGER

01 SKINIRON Device

o Device Main Button



02 3 SAFETY SYSTEM

Safety mechanism for user



- o MICRO PROCESS
Due to Micro process system product will turn off within 5mins when not in use.
- o THERMISTOR
Thermistor maintained at 40-50 degrees
- o BIMETAL
When the temperature reaches over 55 degrees the power will automatically turn-off

! PRODUCT DIRECTION

SKIN MANAGER

01 USE

1. Make sure the product is fully charged.
2. Massage face with cream before use.
3. Power will turn off when adaptor is connected.



Turn on Power.
After charge:
After charging, press and hold the power button to perform both facial and body modes



NUTRITION & ELASTICITY STAGE
[*Press the bottom for a few seconds to start 1st level]
1 TIME PRESS: Face level 1 (Weak:Blue)
2 TIME PRESS: Face level 2 (Strong:Red)
3 TIME PRESS: Body level 1 (Weak:Blue)
4 TIME PRESS: Body level 2 (Strong:Red)
PRESS FOR LONG TIME: Power Off



After Pressing the power, press again for heat massage.
Heat massage helps to open pores to enhance function and body mode helps to enhance body elasticity.
First time: warm 40 degrees Celsius
Press twice: warm 45 degrees Celsius
Press three times: warm 50 degrees Celsius
Thermal function may differ from person to person. check with the palm if your hand for finger before use.

- ▶ According to different skin types repeat process at least two times.
- ▶ Product vibrates when turned on. LED mode automatically turns on.
- ▶ When not in use, product will turn-off within 5mins.

! PRODUCT DIRECTION

SKIN MANAGER

02 FACIAL MASSAGE ORDER



1. Turn on power button



2. Before starting apply nutritional Cream on face and apply Thoroughly and use recommended. After 5mins Facial mode power Will automatically turn off.



3. For stronger massage press one more time (1level:Blue, 2level:Red)




4. Press for a few seconds to turn-off.


03 FACIAL MASSAGE POINTS



Ear: Behind the ear
Before massaging press up and Down behind the ear to fully.

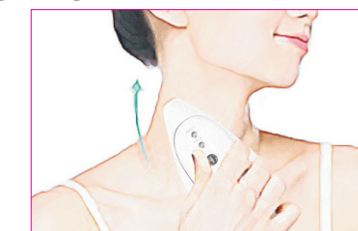

Neck
Use up and down motions in massage motion.

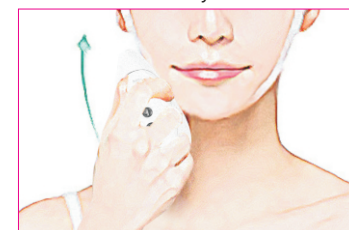

Face Contour
Place on the bottom of the face and move up and down like a real iron.



Forehead
Move back and forth in massage motion

※ Repeat process for the other side.

04 BODY MASSAGE STEPS



1. Press to Body Mode



2. Apply body cream before use and use heat mode, recommended. Gently massage with thermal mode for 5 minutes, the device will turn off.



3. For Stronger massage press one more time
(1 level: Blue, 2 level: Red)



4. Press for a few seconds to turn-off.

05 BODY MASSAGE POINTS



Armpit
From the shoulder to the armpit use sweeping motion.



Arm
Start from the inner side towards the outer elbow area.



Belly
From the upper belly to the bottom use sweeping motions downwards.



Waist
From the top side of the belly use sweeping motions downwards.



Thigh
Use the device in a massage motion similar to when using an iron to spread treatment.



Calf
Use Circular motions up and down.

※ Repeat process for the other side.

01 ADAPTER

◦CHARGER



◦HOW TO CHARGE



#Charge before use.

◦DURING CHARGE

Charging will start with red light,once charge is complete light will change to blue.
(disconnecting the power jack will discontinue the charge)
Charging via usb port may take longer to fully charge.

◦LOW BATTERY

Red light will start to blink

◦ON FULL CHARGE

Product will be available for use for 80 mins. (Lithium battery 3.6v / 2.600mAh)
Please charge fully before use.

※Battery life may vary according to conditions and use.

Adapter Input : C100V-240V 50/60Hz

Output:DC 5V =1.0A

Adapter may vary by country.

01 SKINIRON CAUTION IN USE

- 1.Please read the direction thoroughly before use.
- 2.Use after fully charging.
- 3.Use the proper adapter.
- 4.Remove all jewelry before use.
- 5.Do not use on scars or open wounds.
- 6.Do not use on skins prone to blemishes,allergy,and other related skin problems.Seek medical advice from certified doctors before use, Should skin redden,swell,or itch stop use immediately and seek medical advice from a certified doctor.
- 7.Do not use in wet environment or during shower or bath.
- 8.Take care not to burn skin.
- 9.Do not throw or drop from high places.
- 10.Should there be any sign of defect or doubts in functions please Refer to our sales representatives.
- 11.Product is not for pregnant women,use by child and weak should be supervised.

02 PRODUCT STORAGE & CAUTION

- 1.Please turn off product when not in use.
- 2.Do not disassemble or take apart.
- 3.Do not place in high temperatures or source of heat.
- 4.Product is not water proof,wash with cleaning towel or cloth and do not place under water.
- 5.After fully charging please store with power jack removed.
- 6.Please cleanse head part after use.