# ULTI

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## 1. Product overview

Welcome to Use JULTI beauty instrument, this beauty instrument has the effect of tightening the face and delaying aging! Please read this manual carefully before use and strictly follow the instructions in the manual to operate the beauty instrument

## 2. Technology principle

By transferring 60-65 °C energy waves to the epidermis, dermis and fascia layers of the skin, the skin tissue can produce fine and independent heat condensation points, activate the self-rebuilding ability of the skin, promote the regeneration and recombination of skin collagen and elastic fiber, and achieve the effect of layer by layer tightening and pulling.

### 3.Parameter and list

#### Instrument Parameter



Product name	JULTI	Model	H-033
Power	42W	Input	100V-240V
Voltage	28V	Packing size	265*220*53mm

#### Packing list

1 Host	1 Power adapter	1 Power cable
Cartridges : 4 pcs (Big : 3.0/1.5mm), small :3.0/1.5mm)		
1 pcs Small operation head	1 Manual	Certificate

## 4.Product structure
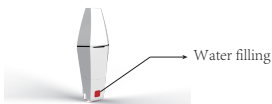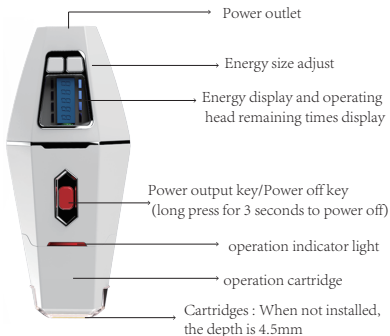
## 5. Operation Steps

1. Clean the operating parts before use.

2. Connect the power cord and press the red button to turn on,

long press for 3 seconds to turn off.

3. Use a large operating head,

select and install a suitable

depth adjustment piece

4. Adjust the energy output, the lowest is 1 and the highest is

8. Press the plus and minus button

to set the increase

or decrease of energy. Adjust the energy according your

feeling (it is recommended to start adjusting from energy 1.)

(1) Energy 1 adjustment, the first row of energy indicator on the

display screen

flashes; Energy 2 adjustment, the first row

of energy indicator on the energy screen

.)

(2) Energy 3 adjustment, the second row of energy indicator on

the display

flashes; Energy 4 adjustment, the second row

of energy indicator on the energy screen is always on.

(3) Energy 5 adjustment, the third row of energy indicator on

the display screen  flashes; Energy 6 adjustment, the third row of energy indicator on the energy screen is always on.

(4) Energy 7 adjustment, the fourth row of energy indicator on

the display  flashes; Energy 6 adjustment, the fourth row of energy indicator on the energy screen is always on

5. Press the energy output button  and the machine starts to

work, and the flashing indicator light indicates  the

output status. Once the energy is outputted,  the

indicator light stops flashing.)

6. After the operation is completed, apply a moisturizing mask on the cleansed face for 15-20 minutes.

**Notice :**

1. The probe must be close to the skin during operation.
2. Forehead, cheekbones, eyes and other areas with less fat, it is recommended to lower the energy.
3. After each operation, the probe needs to be wiped clean with water or saline.

## 6. Operation Method

Medium: gel, aloe vera gel

### Operation method of face:

1. When use big handle, First operate the 4.5mm, then 3.0mm probe, then operate

the 1.5mm probe. (therefore, the operation of probe is from deep to shallow) Point operation according to the direction indicated by the line

2. Operate on face, According to “井” type.

Start from the cheek and divide it into three lines from the inside to the outside for fixed-point operation.

Operate from the chin to the earlobe

Operate from the corner of the mouth to the direction of the ear

Operate from the side of the nose to the temple;

Then from the cheek, divide three lines from bottom to top for fixed-point operation

Operate from the mandible to the inner side of the cheekbone,

Operate from the middle of the mandible to the middle of the cheekbone,

Operate from the masseter muscle to below the temple

## Operation Method for forehead:

1. When use big handle, First operate the 3.0mm probe, then operate the 1.5mm probe. (therefore, the operation of probe is from deep to shallow) Point operation according to the direction indicated by the line

2. Operate from the center of the two eyebrows from up to the fixed point next to the hairline, and then from the upper part of the eyebrows, the peak of the eyebrow, and the tail of the eyebrow to the direction of the hairline.

## Operation Method for neck:

When use big handle, First operate the 3.0mm probe, then operate the 1.5mm probe. (therefore, the operation of probe is from deep to shallow) Point operation according to the direction indicated by the line

2. Move from the top of the collarbone up to the skin, avoiding the position of the Adam's apple;

## 7. Contraindications

People with the following symptoms should be cautious when using this instrument. Please consult a doctor or professional before using it. The details are as follows:

1. For those who have a prosthesis at the operating site, it is recommended to avoid this area when operating the instrument;
2. It is not recommended to use in women's menstrual period;
3. If there is hyaluronic acid or collagen injection in the operation area in the last six months, please let us know;
4. It is not recommended to operate with a serious heart disease and a heart rhythm adjuster;
5. Avoid places such as wounds or purulent acne during operation;
6. It is not recommended to use during skin inflammation or skin allergies;
7. Although it has not been proved that the instrument is harmful to pregnant women, it is not recommended for pregnant women to use it;
8. The instrument cannot be operated on the eyes, nose, ears, mouth, upper eyelids, apple throat, etc.

**Warning: It is forbidden to use it in the bath or showe**

## 8. Precautions after operation

1. Do not wash your face with hot water within three days (you can use warm or cold water to wash your face);
2. Do not soak in hot springs or go in or out of high-temperature places within three days;
3. Strengthen the moisturizing and sun protection of the operating parts;
4. It is recommended to apply a moisturizing mask every day within three days of use;
5. Eat less spicy food within three days, smoking, drinking, staying up late is not recommended

## 9. Troubleshooting

Fault type	Troubleshooting method
After plugging in the power, the instrument does not respond	First check whether the power strip is powered on, and replace the power strip or socket

## 10. Instrument maintenance

1. The instrument must use a three-pronged plug with a grounding pin, and ensure that the power socket of the instrument is well grounded;
2. Ensure that the power supply is stable and adaptable. If the local power supply voltage is unstable, it is recommended that the user add a regulated power supply with matching power;
3. The instrument should not be placed in a humid place or near a water source, and the instrument should not be directly exposed to sunlight;
4. Keep the instrument away from strong pyrogens or strong radiation sources (disinfection cabinets, induction cookers, microwave ovens, etc.), which will affect the normal use of the instrument;
5. The instrument must use the original qualified accessories provided by the manufacturer;
6. After each operation, the operating handle must be cleaned and disinfected and kept properly;
7. When the instrument is not in use, please turn off the power, unplug the power cord and properly place various accessories of the instrument.

# 11. Warranty card

To users:

Thank you for your support and trust!

xx If any damage caused by the manufacturing process or components occurs within one year of the date of purchase, a free warranty is promised (note: man-made damage must be repaired by individuals)

Warranty

1. Warranty date within one year, starting from the date of purchase
2. Please present the proof of purchase or purchase record for the warranty product;
3. After the warranty date, repair and replacement parts will be charged at the price of replacement parts;
4. In one of the following situations, free warranty is not implemented:
  - 1) The handle and body were dismantled during use, causing the instrument not to be used normally;
  - 2) Unable to provide the following voucher or purchase record

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