


Smil body User manual



Contents



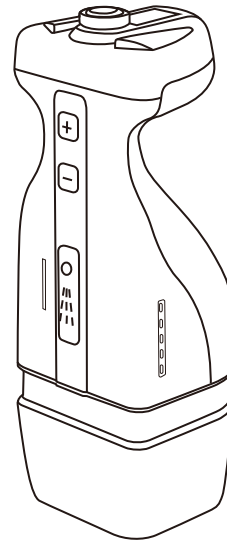


Product introduction	1
Product list	2
Product parameter	2
Instrument instruction	3
Steps	5
Manipulation.....	7
Operation contraindications.....	8
Instrument maintenance.....	8
Common faults and troubleshooting.....	9
Warranty card.....	9

Product introduction

The rotary fat-reducing mechanical wave is used to operate at 8-13MM deep into the skin, which interferes with the hot cells of the subcutaneous fat tissue, uses the heat to permanently destroy the subcutaneous fat, directly reduces the excess fat cells, and stimulates the contraction and increase of collagen. Thick, stimulate elastic fibers and destroy fat cells, while the surrounding important structures such as blood vessels and connected tissues are still intact.



Product List


Product name	Host	Power cable	Sticker	Instruction	Warranty card
Product Image					
Quantity	1	1	2	1	1

Product parameter

Product: Smil body	Depth:8mm~13mm
Model no: H-028S	Machine size: 68*68*198mm
Power: AC110V~240V/50Hz~60Hz	Package dimension: 525*343*432mm
Output power: 24W	weight: 1.0kg
Frequency: 4MHz	

Instrument instruction

(Product top view)

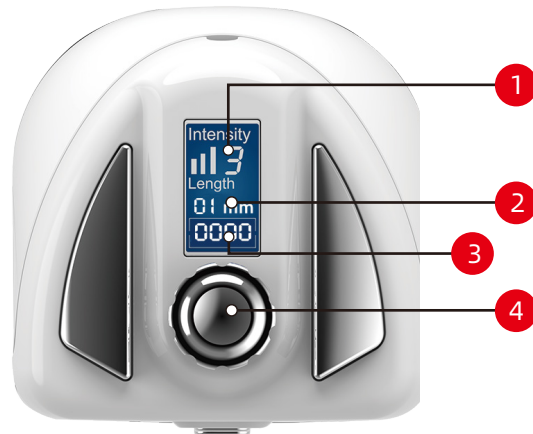

- (1) The current energy intensity is displayed during operation. The maximum energy is 9 gears and the lowest is 1 gear. It can be adjusted according to the feeling during operation.
- (2) The distance between the energy output point and the point does not need to be adjusted
- (3) It will automatically record every operation, no need to adjust
- (4) Depth control valve, respectively 8MM and 13MM depth, fat thickness below 20MM use 8MM depth operation, fat thickness over 20MM use 13MM depth operation




- (5) During operation, the energy intensity can be adjusted through the +--button
- (6) Depth scale, the position of the small red line aligned with 8 represents the depth of 8MM, and the alignment 13 represents the depth of 13MM, which can be adjusted by rotating the (4) depth adjustment valve
- (7) Operation display
- (8) Power cord jack
- (9) ON/OFF
- (10) Head protective cover, after the operation is finished, the instrument head should be cleaned and used to avoid the abrasion and scratches caused by the head touching other objects.

Steps


1. Clean the operating part and apply the medium firming emulsion. The medium needs to be coated with a thickness of about 3mm, and it cannot be coated too thin. The head should also be coated with the medium, adjust the depth and energy, and start the operation.



2. During the operation, the head must be close to the skin. Do not operate in the air. Do not move the operating head during the energy output. After one shot is finished, move the distance 3-4 mm and operate the next one close to the edge of the previous shot. hair.

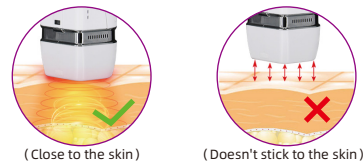

3. Two shots cannot be operated continuously at the same position, and the operation can be repeated 3-5 times back and forth.

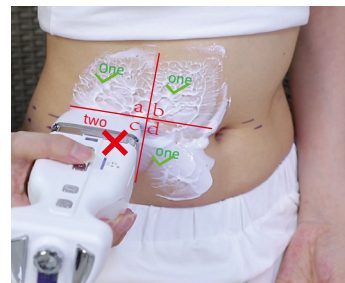

- The operation time of one part is 45-60 minutes
- At the end of the operation, massage by hand to absorb the firming cream, then wipe it clean with a wet wipe

Attention:

- Because the body and facial fat distribution are different, the instrument cannot operate the face







- The delivered compacting cream should be applied to the operation part or wrist for allergy test before operation. After 3 minutes of massage, wipe it off. After 5 minutes, there is no discomfort such as itching, red bumps and other allergic symptoms in the test area, you can use this product to operate. If you have symptoms such as itching, you need to replace other media. Note: It is normal for the firming cream to have a slight hot sensation after operation. Pay attention to the difference between mild hot and itchy.

Course of treatment:

- 10~15 times are recommended as a course of treatment. Operate the first 3 times a day, and operate it once in 3-5 days

Manipulation

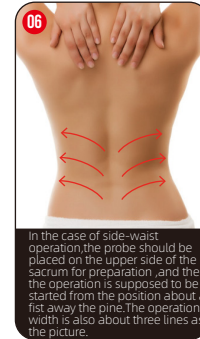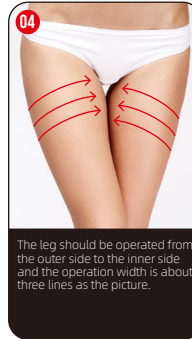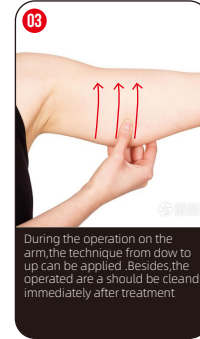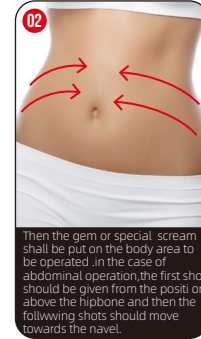

■ Precautions:

1. Do not take a bath within two hours after operation
2. Do not take a hot bath within 24 hours, use warm water not exceeding 40 degrees
3. After the operation, keep drinking more than 2000ML of water every day, 4 bottles of regular 550-600ML mineral water
4. Light diet, eat less sweets, and not overeating
5. Regular work and rest, avoid drinking, staying up late, supper and other habits during operation
6. The instrument needs to add water once every 4-5 times of use, one part is added once, and pure water or mineral water is added, not tap water.

Operation contraindications

- (1). It is recommended to avoid areas with wounds or scars on the skin;
- (2). Use hormonal drugs with caution;
- (3). Prohibited for pregnant women, lactation and menstrual periods;
- (4). Built-in metal pieces in the body (such as pacemakers, etc.)
- (5). People who are extremely weak
- (6). Within three months after the female wears the birth control ring
- (8). Females cannot operate their abdomen within three months after normal delivery, and do not operate their abdomen within half a year after caesarean section
- (9). Skin during allergy
- (10). Severe diseases, such as heart disease, gastric ulcer, severe gastric disease, duodenal ulcer, high blood pressure, diabetes, etc.
- (11). Medical plastic parts or parts and artificially filled substances need to consult a doctor first.



Warning sign Note: Do not use while bathing or showering.

Instrument maintenance

1. The instrument must use a three-pronged plug with a grounding pin, and ensure that the power socket of the instrument is well grounded;
2. Ensure that the power supply is stable and adapted. If the local power supply voltage is unstable, it is recommended that the user add a regulated power supply with matching power;
3. The instrument should not be placed in a humid place or water source accessories, and the instrument should not be directly exposed to sunlight;
4. The instrument should not be too close to a strong heat source or a strong radiation source (sterilization cabinet, induction cooker, microwave oven, etc.), which will affect the normal use of the instrument;
5. The instrument must use the original qualified accessories provided by the manufacturer;
6. It is recommended to use the treatment supporting products provided by the manufacturer. If you use outsourcing products provided by non-manufacturers, you must first consult the manufacturer and obtain the manufacturer's consent before use. For outsourcing products, it is recommended to use related products from formal channels;
7. After each operation, the function handle and function head must be cleaned and disinfected and kept properly;
8. When the instrument is not in use, turn off the power, unplug the instrument, and place various accessories of the instrument. If possible, cover the instrument with a dust cover.

Common faults and troubleshooting

Faults	Trouble shooting
No react when turn on power	To check if there is electricity in socket, change other socket to have a try.
Head's water leaking	Stop using it, change the plastic cover or return back the machine to the manufacturer.
Please note: This machine is precision instrument and it is forbidden to be used in the humid environment such as bathroom. Please stop using it and return back to the manufacture once it gets into water.	

Warranty card

To users:

Thank you for your support and trust!

If Smil body is damaged due to manufacturing process or components within one year from the date of purchase, a free warranty is promised (note: man-made damage must be repaired individually)

Warranty

-The warranty date within one year starts from the date of purchase;

-Please show proof of purchase or purchase record for warranty products;

-After the warranty date, repair and replacement parts will be charged at the price of replacement parts;

-In one of the following situations, free warranty is not implemented:

-The handle and body were dismantled during use, causing the instrument not to be used normally;

-Unable to provide the following proof or purchase record.

Model No.:

Purchase date:

Client:

Phone No.:

Fault description