



Guangzhou Huaxi C&T Co.,Ltd.

Add: Room 1502 Hongmian Business Center, No.10 Hualei Road, Liwan District, Guangzhou, Guangdong Province, China

E-mail: info@choicybeauty.com Tel:+86-20-38263199 Website: www.choicybeauty.com

Choicy Beauty Equipment Series

-- DPL Multi Service Platform



Model: CB-1201A

Product Description




System Characteristics

1. Machine is used the multiple spectrum and narrow spectrum design for a variety of facial physical problems and effective simulation of laser efficacy, so you can enjoy safety and effectiveness unity.
2. More lasting reliable quality, the main parts are used the same as the US and Israel machines' using.
3. More efficient cooling: the safer treatment can reduced guests pain significantly
4. More high security: multiple fault diagnosis can help the body multiple self-protection.
5. Improve the hand-piece performance: you can use smaller energy parameters and fewer number of treatments to achieve a more lasting and significant effect.

Treatment Range

1. Face rejuvenation, improve overall facial / skin texture (550-950nm)
2. Shrink pores (550-950nm)
3. Cure freckles, sunburn, age spots and other pigmentation (DPL 550-650nm)
4. Cure telangiectasia, acne redness, facial flushing, port wine stains and other vascular diseases (DPL 500-600nm)
5. Acne treatment
6. Fast and effective hair removal

Before and After Treatment

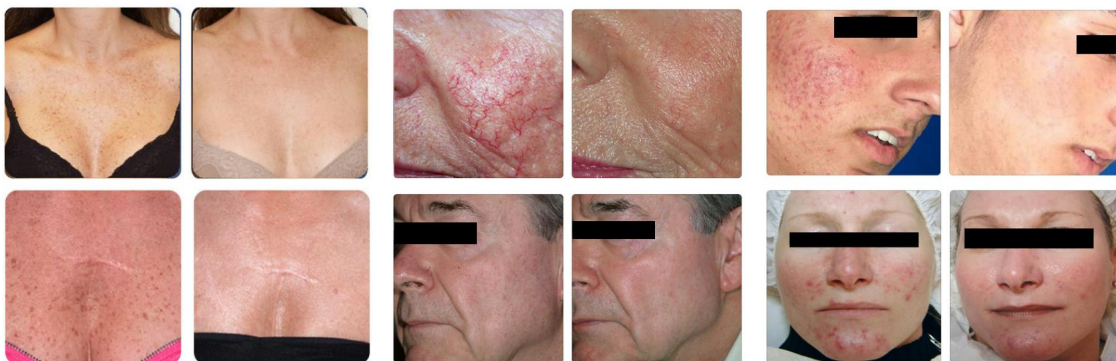



Guangzhou Huaxi C&T Co.,Ltd.

Add: Room 1502 Hongmian Business Center, No.10 Hualei Road, Liwan District, Guangzhou, Guangdong Province, China

E-mail: info@choicybeauty.com Tel:+86-20-38263199 Website: www.choicybeauty.com

Technical Specifications

Standard wavelength	550-950nm; 650-950nm
Optional wavelength	420-950nm; DPL 500-600nm; DPL 550-650nm
Energy density	10-35j/cm ²
Spot size	10 x 35mm
Cooling system	Semiconductor cooling & water cooling & air cooling
Treatment model	
Simple model	
Pulse number	1
Pulse width	10ms, 12ms, 15ms
Pulse delay	Intelligentize delay
Professional model	
Pulse number	1-4
Pulse width	10ms
Pulse delay	100ms

NOTE:

The detailed instruction will be sent with the equipment.